





Meditech Kft.
 Budapest 1184 Mikszáth Kálmán u. 24.
 Tel. 36(1)280-82-32
 Fax. 36(1)282-93-88

Patient data			
name	Demo card(X)plore		
ID	0001-0006	date of birth	1961. 10. 28.

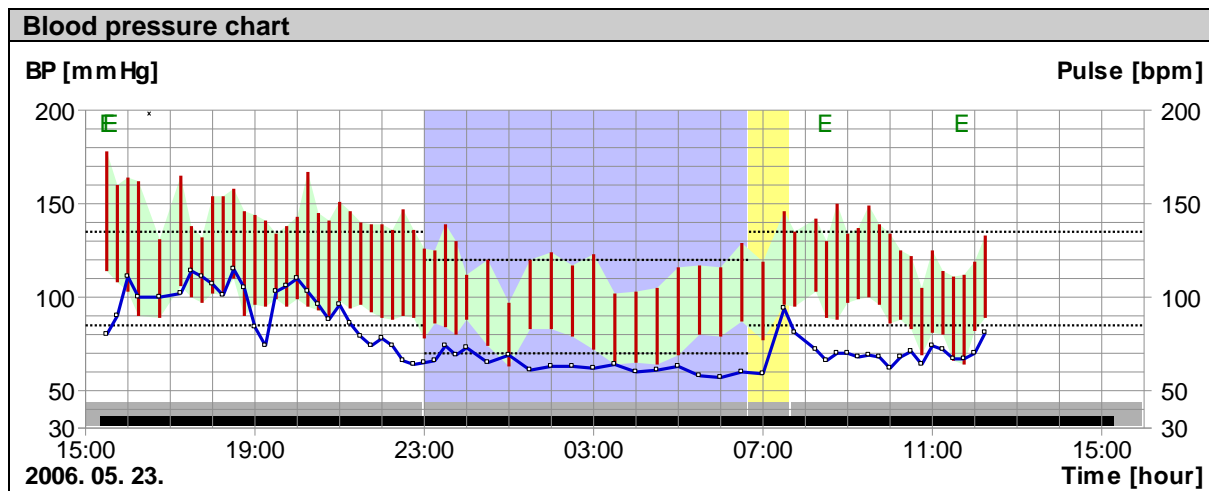


ABPM basic data (2006. 05. 23. 15:20 - 2006. 05. 24. 15:20)				
period	Overall	Morning	Day	Night
time range	Overall	06:40 - 07:40	07:40 - 23:00	23:00 - 06:40
SBP/DBP maximum limit value	~127/77~	135/85	135/85	120/70
SBP/DBP weighted average	131/86	131/86	140/93	116/76
Pulse weighted average	76	73	84	63
readings	66	4	47	18
success ratio %	96	100	94	100
SBP/DBP D/N index %	17/18			
morning surge	21			

Hypertension summary (2006. 05. 23. 15:20 - 2006. 05. 24. 15:20)				
period	Overall	Morning	Day	Night
time range	Overall	06:40 - 07:40	07:40 - 23:00	23:00 - 06:40
SBP/DBP maximum limit value	~127/77~	135/85	135/85	120/70
SBP/DBP weighted average	131/86	131/86	140/93	116/76
MAP weighted average	101	101	108	89
PP weighted average	45	45	48	41
Pulse weighted average	76	73	84	63
SBP/DBP PTE %	50/78	38/50	62/85	30/67
SBP/DBP load mmHg	148/199	99/113	206/217	52/172
readings	66	4	47	18
SBP high readings	38	1	30	7
DBP high readings	53	3	39	13
SBP/DBP D/N index %	17/18			
SBP/DBP D/N mean difference	24/17			
morning surge	21			


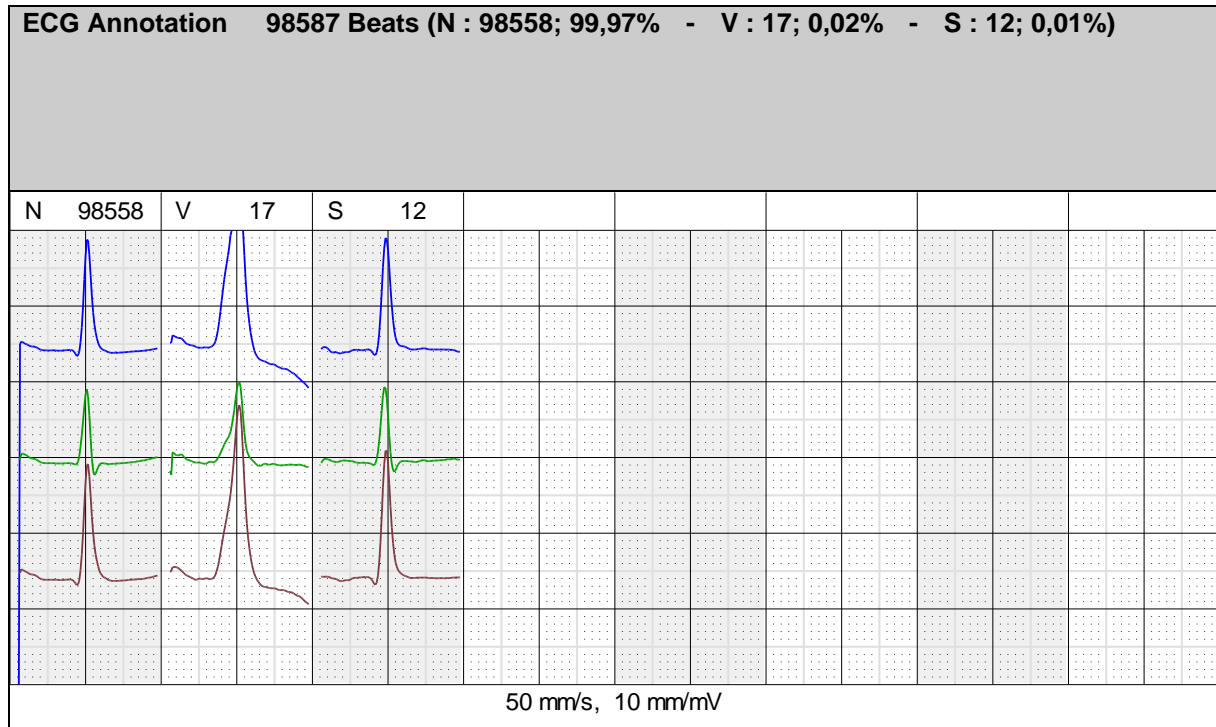


Holter summary					
The measurement started	2006. 05. 23. 15:25	Total beats	98587		
Total measurement time	20 h 51 min	Normal	98558	(100,0%)	
Valuable length:	20 h 50 min	VES	17	(0,0%)	
		SVES	12	(0,0%)	
Channels / sampling rate	3 / 300 Hz	PM	0	(0,0%)	
Total lead off / noise total	00 h 00 min	Other	0	(0,0%)	
Pacemaker spike detection	disabled	Artefact	0	(0,0%)	
HR and RR intervals	Average	Max	Min	Max time	Min time
HR [1/min]	78	130	55	2006. 05. 24. 07:24	2006. 05. 24. 06:01
RR filtered [ms]	760	1293	446	2006. 05. 24. 06:09	2006. 05. 24. 07:24
RR unfiltered [ms]	760	1293	423	2006. 05. 24. 06:09	2006. 05. 24. 12:13
HRVti: 20 based on this the risk factor is moderate					
SDNN: 157 based on this the risk factor is low					
Tachy / brady / pause	All events		Longest event		
	Number	Total length	Time	Length	
Tachycardia	97	4 h 3 min 48 sec	2006. 05. 23. 17:23	45 min 36 sec	
Bradycardia	0	0 sec			
Pause (2.5 sec)	0				
Rhythm disturbances					
Supraventricular ectopy:	rarely occurring SVES beats				
Ventricular ectopy:	rarely occurring VES beats and V couplets				
ventricular fibrillation, atrial fibrillation not found					
Other data					
QT maximum:	430 ms (06:33)	QTc maximum:	499 ms (17:50)		
Total ischaemic burden:	stA=00:02	stB=00:00	stC=00:00		
Blood pressure data					
Average blood pressure [mmHg]		Total 24 h:	131/86	Day:	140/93
				Night:	116/76
PTE (Hypertensive time index):	40/58 %		Diurnal index(Day/Night index): 17/18		

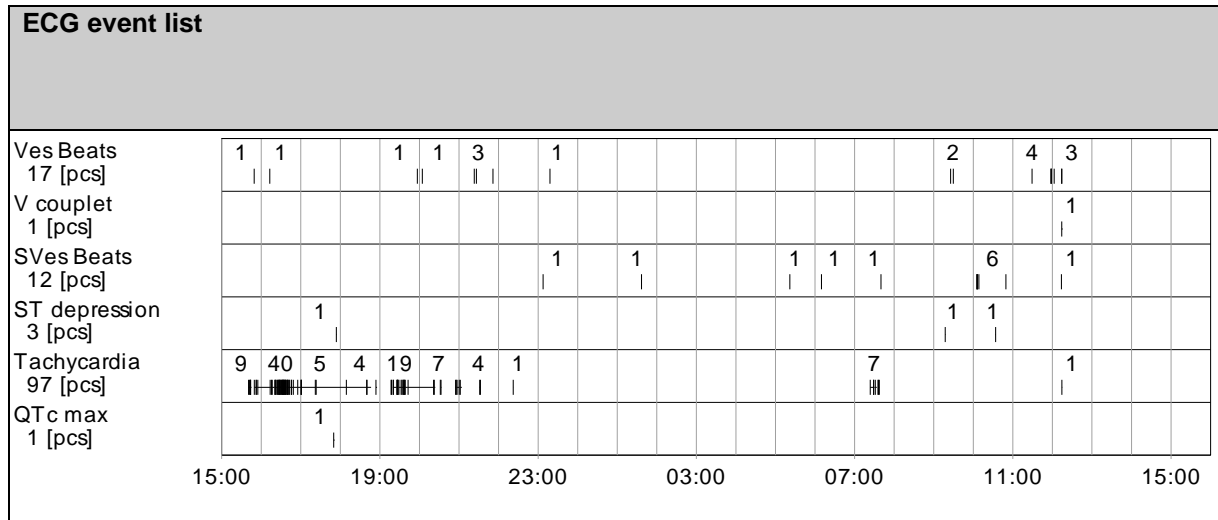

HRV statistics

Risk 24 hour SDNN: 157	HRV statistics (Custom) 2006. 05. 23. 15:20 - 2006. 05. 24. 15:20					
low (above 100)	SDANN	149 ms	SDSD	11 ms	RR count	98580
	SDNN index	34 ms	NN50	139 pc	RR average	759 ms
mid	NN count	98555	pNN50	0 %	RR std.dev. (SDRR,CLV)	157 ms
	NN average	759 ms	HRVti	20	rMSSD	12 ms
high (below 50)	NN std.dev. (SDNN)	157 ms	TINN	314 ms		

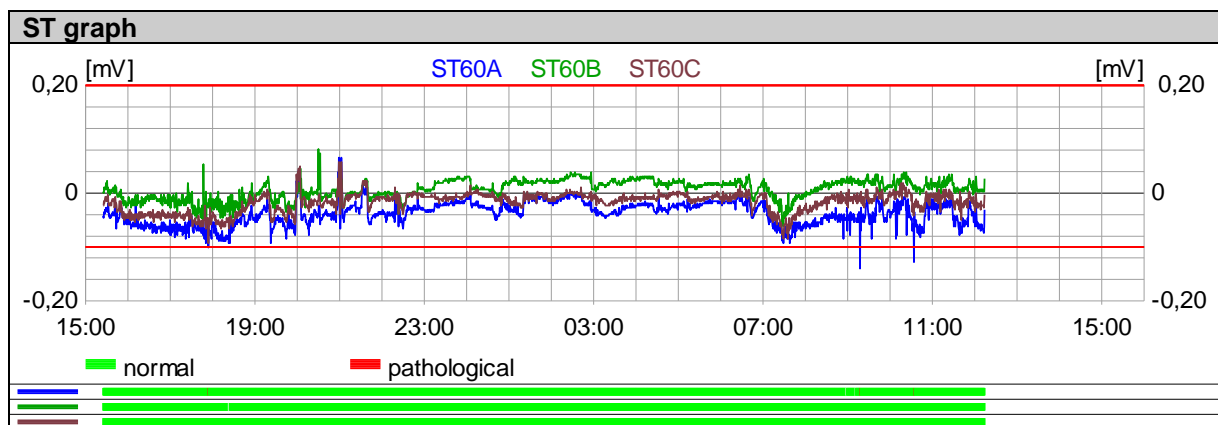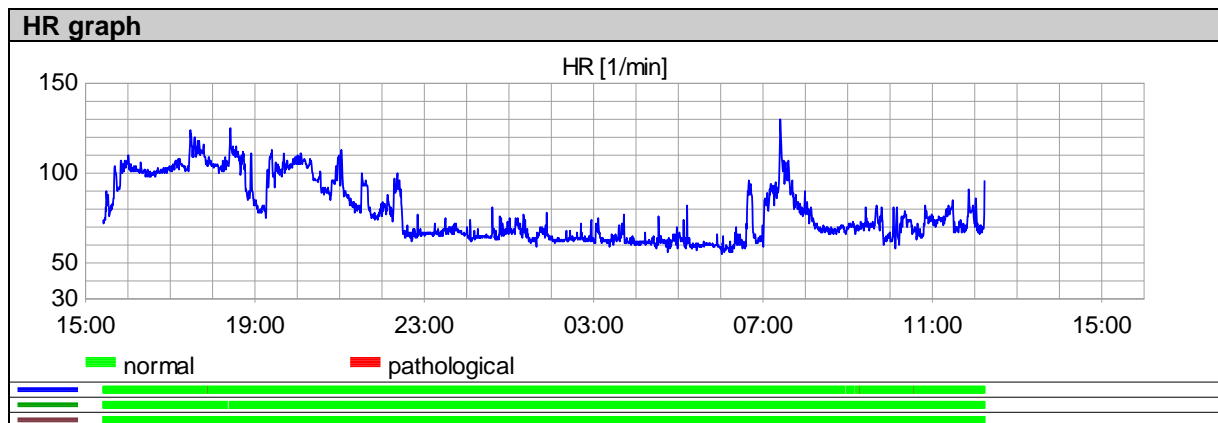

ST event list								
start	ch	duration [min]	STmax [mV]	TIA [mV*min]	HRstart [1/min]	HRmax [1/min]	HR [%]	Type
2006. 05. 23. 17:54	A	0,5	-0,10	-0,08	107	109	1	depression
2006. 05. 24. 09:17	A	0,5	-0,14	-0,09	69	70	1	depression
2006. 05. 24. 10:34	A	0,5	-0,13	-0,08	68	68	0	depression

