




## All-day ECG Monitors

Hearts Deserve  
the Best Holters

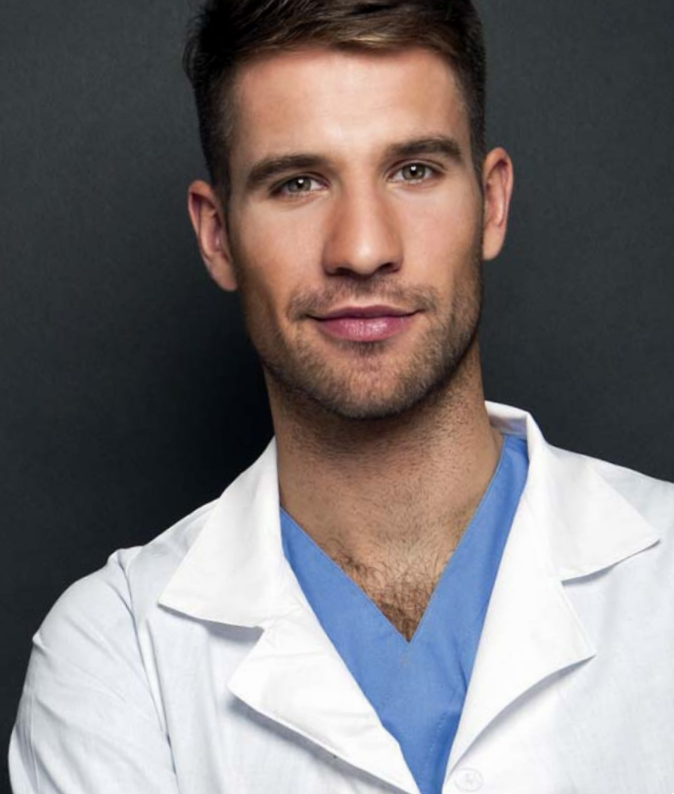

# CardiUP!12

3/12-channel Holter ECG



## CardiUP! benefits

- stay connected up to 7 days for better cardiovascular analysis
- 3/12 or 3 recording channels
- large LCD screen for real time ECG display
- event button
- both manual & software programming
- mini USB connection
- new, sensitive PM detection
- multifunctional carrying case
- free, all-inclusive software

## CardiUP!3

### 3-channel Holter ECG




# CardioMera

basic, multichannel Holter ECG

## CardioMera benefits

- 1-, 2-, 3- or 5-channel recording
- individual, easily replaceable patient cables
- small size & light weight
- event button
- free, all-inclusive software

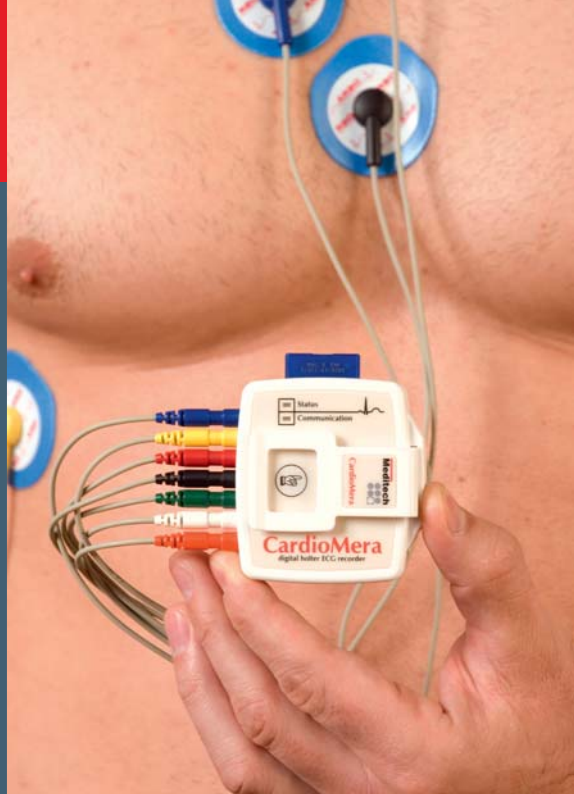

<b>Specifications</b>	<b>CardiUP!12</b>	<b>CardiUP!3</b>	<b>CardioMera</b>
dimensions	67*93*21 mm		66*59*17 mm
weight (excl. batteries)	92g		50g
memory	internal		external (a card reader unit is needed for retrieving data)
data transfer	USB 2.0 (mini USB)		USB optical cable
power supply	1 AA battery		1 AAA battery
LCD display	48*28 mm		–
number of keys	4 (event, up, down, select)		1 (event)
number of channels	3/12 independent channels (7 or 10 leads)	3 independent channels (7 leads)	1, 2, 3, 5 channels (7 leads)
operation	up to 168 hours of continuous ECG		up to 96 hours of continuous ECG
sampling, A/D resolution	300, 600 or 1200 Hz, 12 bit		600 or 1200 Hz, 12 bit
storage quality	300, 600 or 1200 Hz		150, 300 or 600 Hz
conformity certificates	CE		

## Holter ECG REPORT



### PATIENT DATA

name	0001_0001	Meditech	2001.01.01.
ID		date of birth	

### HOLTER SUMMARY

The measurement started	2018.03.09. 11:57	Total beats	96311	
Total measurement time	24 h 00 min	Normal	10553	(11,0%)
Valuable length:	24 h 00 min	VES	127	(0,1%)
		SVES	223	(0,2%)
Channels / sampling rate	12 / 300 Hz	PM	85408	(88,7%)
Total lead off / noise total	00 h 00 min	Other	0	(0,0%)
Pacemaker spike detection	enabled	Artifact	0	(0,0%)


16:46:48

A+  
HR=82

HR

| 74 | 76 | 74 | 104 | 127 | 82 |

### HRV TACHOGRAM

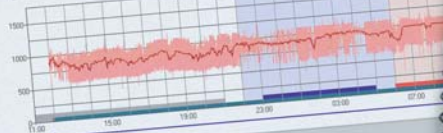

### HRV BAR CHART



## CardioVisions free software for ECG analysis

Specifications	CardiUP!2	CardiUP!3	CardioMera
accuracy	99,90%		
beat-by-beat ECG annotation		✓	
event distribution		✓	
ST analysis		✓	
QT analysis		✓	
PQ analysis		✓	
PM analysis		✓	
time domain HRV analysis		✓	
frequency domain HRV analysis		✓	
arrhythmia analysis		✓	
graphical and tabular display of data		✓	
customizable screens		✓	
customizable & printable report		✓	
personal & network use		✓	



Meditech Ltd • H-1184 Budapest, Mikszáth Kálmán u. 24.  
phone: +36 1 280 8232, +36 1 280 8233 • email: [meditech@meditech.eu](mailto:meditech@meditech.eu)  
[www.meditech.eu](http://www.meditech.eu)